

INDUSTRIAL SOLUTIONS  
**A.P.I.P.S**

AEROSPACE / INDUSTRIAL  
BIOPHARMACEUTICAL  
BIOTECHNOLOGY  
MEDICAL DEVICE  
PHARMACEUTICAL

# IPA 7030 STERISPRAY®

70% USP Isopropyl Alcohol, USP  
30% Water for Injection, USP by Volume



Aerospace / Life Sciences / Microelectronics



## IPA 7030 STERISPRAY®


70% USP Isopropyl Alcohol, USP  
30% Water for Injection, USP by Volume

STERIL, VALIDATED AND CERTIFIED

# Features

STERILE, READY-TO-USE 70% ISOPROPYL ALCOHOL (IPA) SOLUTION  
DESIGNED ESPECIALLY FOR CLEANROOM USE

Made with Water-For-Injection (WFI), filtered to 0.2  $\mu\text{m}$ , individually double-bagged and gamma-irradiated to ensure sterility, It can be used to sanitize gloves, clean surfaces and equipment found in cleanroom areas within pharmaceutical, biotechnology, and medical device facilities or anywhere the use of a sterile, non-pyrogenic alcohol solution is required.

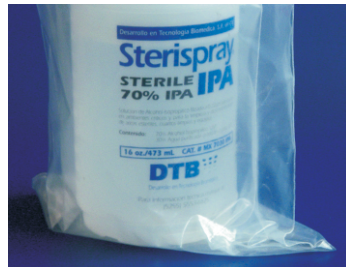

## IPA 7030 STERISPRAY®

70% USP Isopropyl Alcohol, USP  
30% Water for Injection, USP by Volume

STERIL, VALIDATED AND CERTIFIED

# Features

IPA 7030 STERISPRAY®, Made according to USP specifications and packaged in a Class 100 cleanroom. Each lot undergoes a bacterial endotoxin test to ensure that the product meets or exceeds "Water-for-Injection" quality in endotoxin limits and each lot also undergoes a USP 14-day sterility test.

Each case is shipped with lot specific documentation which details QC, irradiation, sterility and bacterial endotoxin limits. Available in 16 oz and 32 oz non-aerosol trigger spray bottle.

ALL LOTS ARE GUARANTEED NOT TO EXCEED THE ENDOTOXIN LEVEL CRITERION FOR WFI .

- Filtered to 0.2 microns.
- Gamma-irradiated to a sterility assurance level (SAL) of  $10^6$
- Each lot is tested for sterility via USP 14-day sterility test.
- Each lot is tested for bacterial endotoxin via the LAL test prior to release
- Bottled and packaged in a Class 100 cleanroom.
- Individually double bagged in a solvent resistant bag.
- Assembled and ready to use.
- A copy of lot specific documentation, which details QC, irradiation, sterility, and bacterial endotoxin limits, is attached to every case.



## IPA 7030 STERISPRAY®

70% USP Isopropyl Alcohol, USP  
30% Water for Injection, USP by Volume

STERIL, VALIDATED AND CERTIFIED

## Benefits

- Ready to transfer and use in critical areas.
- Lot specific documentation, which details QC, irradiation, sterility, and bacterial endotoxin limits.
- Expiration date and lot number on each bottle and each box.
- Quality Control and Validation processes only once.



Graduated Spray bottle for precise dispensing.



## IPA 7030 STERISPRAY®

70% USP Isopropyl Alcohol, USP  
30% Water for Injection, USP by Volume

STERIL, VALIDATED AND CERTIFIED

# Manufacturing vs Purchasing

Minimal manufacturing processes required to obtain this product with a validation certificate involve huge diverse and hidden costs that their allocation dilutes the real financial impact on the company, hence the end user will allocate unreal costs that will not be the real and much higher ones.

Using IPA 7030 STERISPRAY® is easy, convenient and manageable and generates substantial savings.



Adjustable Non-Aerosol Trigger Spray  
to choose a stream for precise application  
or coarse spray for large areas.

## IPA 7030 STERISPRAY®

70% USP Isopropyl Alcohol, USP  
30% Water for Injection, USP by Volume

STERIL, VALIDATED AND CERTIFIED

## Form of use protocol

Use in well-ventilated area. Hold spray bottle six to eight inches from surface and spray area until it is completely covered with solution. Or wipe IPA 7030 STERISPRAY onto surface using cleanroom grade wipers. Allow to air dry.

### BENEFITS OF USING IPA 7030 STERISPRAY

Ready to transfer and use in critical areas.

Lot specific documentation, which details QC, irradiation, sterility, and bacterial endotoxin limits.

Expiration date and lot number on each bottle and each box. Quality Control and Validation processes only once.



Non-Aerosol Fully assembled  
ready-to-use Trigger Spray/Bottle.



## IPA 7030 STERISPRAY®

70% USP Isopropyl Alcohol, USP  
30% Water for Injection, USP by Volume

STERIL, VALIDATED AND CERTIFIED

# Request information

PART NUMBER	DESCRIPTION	PACKAGING
MX7030	IPA 7030 STERISPRAY® 70% Isopropyl Alcohol, USP 30% Water for Injection, USP 16 oz Bottle	12 bottles/case



## IPA 7030 STERISPRAY®

70% USP Isopropyl Alcohol, USP  
30% Water for Injection, USP by Volume

STERIL, VALIDATED AND CERTIFIED

Aerospace / Life Sciences / Microelectronics

# Contact information

Quotes and orders from around the world

APPS INC.  
appsindustrial.com  
+1 (818) 930 0779  
[Toll free]: +1 818-5743-345  
contact@appsindustrial.com

